


HORIBA LA-300 for Windows (TM) Ver. 3.71

LA-300 system for Windows

ファイル名	: 棒焙じ	比表面積	: 2148.3 (cm ² / cm ³)
ID#	: 201103081441063	メジアン径	: 64.9949 (μm)
循環速度	: 2	算術平均径	: 88.2811 (μm)
超音波	: OFF	モード径	: 108.4784 (μm)
透過率 (レーザー)	: 92.3 (%)	残差Rパラメータ	: 1.2318E-1
分布形態	: 標準	バリテーション結果	:
相対屈折率	: 1.30-0.00i		
表示グラフ軸	: LogX-LinY		
物質名	:		
原料	:		
分散処理	:		
備考	: 2011 2/19		

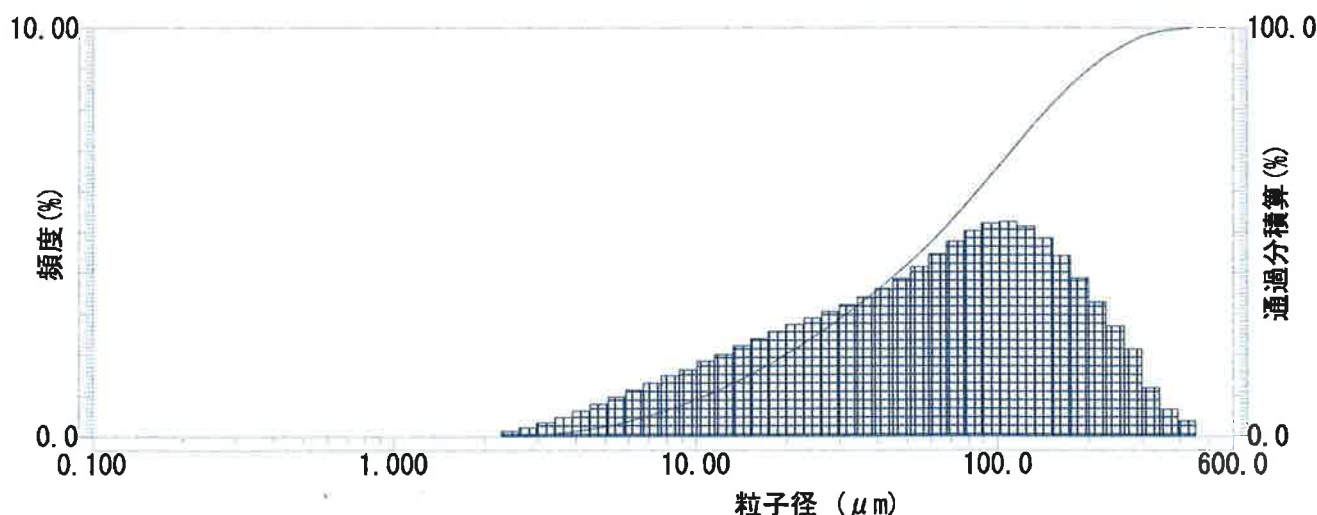

No.	粒子径 (μm)	頻度 (%)	積算 (%)	No.	粒子径 (μm)	頻度 (%)	積算 (%)	No.	粒子径 (μm)	頻度 (%)	積算 (%)
1	0.115	0.000	0.000	30	5.867	0.954	3.501	59	300.518	2.129	97.795
2	0.131	0.000	0.000	31	6.720	1.123	4.624	60	344.206	1.183	98.978
3	0.150	0.000	0.000	32	7.697	1.302	5.926	61	394.244	0.657	99.635
4	0.172	0.000	0.000	33	8.816	1.477	7.403	62	451.556	0.365	100.000
5	0.197	0.000	0.000	34	10.097	1.625	9.029	63	517.200	0.000	100.000
6	0.226	0.000	0.000	35	11.565	1.830	10.858	64	592.387	0.000	100.000
7	0.259	0.000	0.000	36	13.246	1.995	12.854				
8	0.296	0.000	0.000	37	15.172	2.207	15.061				
9	0.339	0.000	0.000	38	17.377	2.371	17.432				
10	0.389	0.000	0.000	39	19.904	2.559	19.991				
11	0.445	0.000	0.000	40	22.797	2.727	22.718				
12	0.510	0.000	0.000	41	26.111	2.878	25.596				
13	0.584	0.000	0.000	42	29.907	3.037	28.633				
14	0.669	0.000	0.000	43	34.255	3.203	31.837				
15	0.766	0.000	0.000	44	39.234	3.389	35.225				
16	0.877	0.000	0.000	45	44.938	3.596	38.822				
17	1.005	0.000	0.000	46	51.471	3.843	42.665				
18	1.151	0.000	0.000	47	58.953	4.135	46.800				
19	1.318	0.000	0.000	48	67.523	4.452	51.252				
20	1.510	0.000	0.000	49	77.339	4.761	56.013				
21	1.729	0.000	0.000	50	88.583	5.024	61.037				
22	1.981	0.000	0.000	51	101.460	5.198	66.235				
23	2.269	0.000	0.000	52	116.210	5.242	71.477				
24	2.599	0.129	0.129	53	133.103	5.122	76.599				
25	2.976	0.214	0.343	54	152.453	4.836	81.435				
26	3.409	0.326	0.669	55	174.616	4.404	85.839				
27	3.905	0.465	1.134	56	200.000	3.866	89.705				
28	4.472	0.623	1.757	57	229.075	3.277	92.982				
29	5.122	0.790	2.547	58	262.376	2.684	95.665				